

APTAMER INFORMATION A3XHX

Description:

- **Identifiers:** A3XHX
- **Number of DNA nucleotides:** 70 bases
- **Target for selection:** Recombinant extracellular domain of CD20

Aptamer was selected from a randomized 32-mer library against the CD20 protein. Proprietary methods were then used to select the aptamer.

Aptamer folding instruction before use:

Once the aptamer is in its working concentration, it needs to be heated to 85-90 °C for 2 minutes, and then cooled to room temperature before use.

1b. Validation data with CD20 by MST (Microscale Thermophoresis) method:

- **Buffer used for validation:** 1X PBS, pH = 7.4, 1 mM MgCl₂
- **Average K_d = 253 ± 24.5 nM**

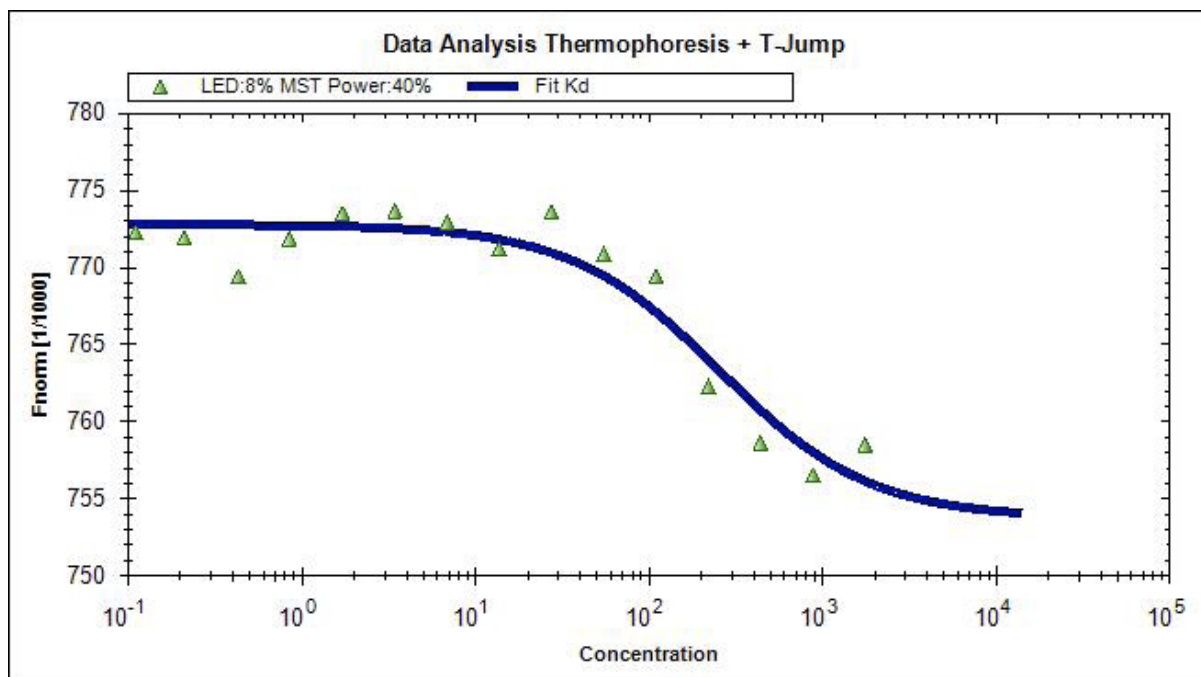


Figure 1. Aptamer-CD20 binding. The thermophoresis is plotted versus the titrated CD20 concentration.